



## ***The triple beam that truly has it all!***


Accurate, dependable and stylish; the OHAUS Triple Pro<sup>®</sup> mechanical balance offers it all: The superior quality and durability you've come to expect from OHAUS, plus all the extras to satisfy every classroom activity requiring a balance. Just like the original Triple Beam balance from OHAUS, the Triple Pro balance gives you the confidence to know you have the accuracy and convenience of a toploading balance while allowing your students to gain hands-on experience through manipulation of the poises, attachment weights and the rod & clamp assembly (all included). The OHAUS Triple Pro balance is the perfect learning tool to help students visualize mass measurement principles and to emphasize a hands-on approach to learning.

### ***Standard Features Include:***

- Unique, easy grip carrying handle is ergonomically designed and perfectly balanced to allow for effortless movement around the classroom
- Two security features are included to prevent theft; a security loop and a Kensington-type lock (cables sold separately)
- Includes a convenient on-board storage compartment for the included attachment weight set
- A built-in transportation/storage lock helps to protect the longevity by minimizing wear on the knives and bearings, and reducing the possibility of bending the damper vane
- Includes rod & clamp assembly for below-balance weighing (specific gravity determination)
- Backed by an incredible 10 year limited warranty that only OHAUS can provide

# Triple Pro<sup>®</sup> Mechanical Balance

## Standard Features:

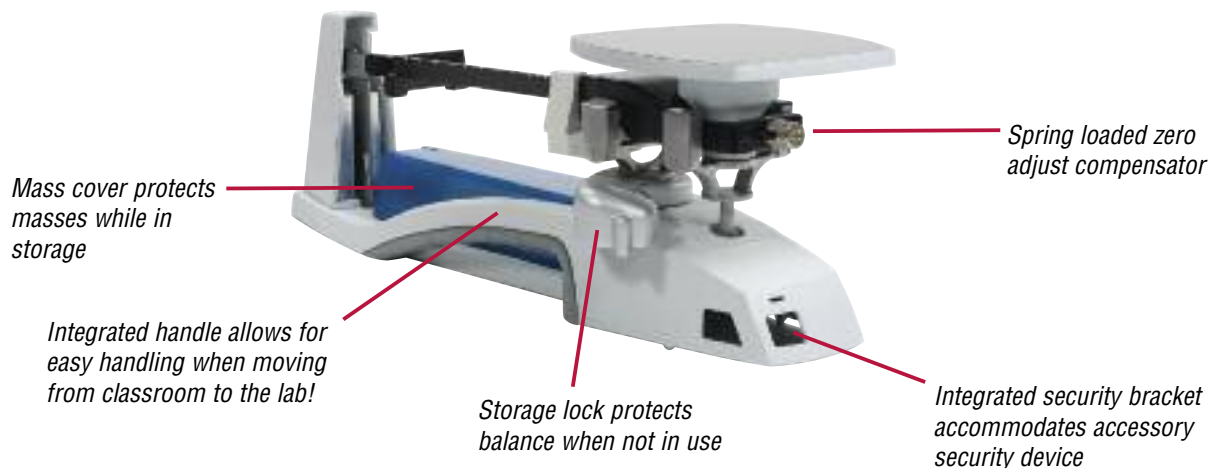

Included Rod & Clamp accessory for specific gravity experiments



Three masses included with balance



Security Device accessory (Lock & cable sold separately)

## Specifications

Model	TP2611
Capacity	2610 g (when using the included attachment weight set)
Readability	0.1g
Calibrations:	
Front beam	10 g x 0.1 g
Middle beam	500 g x 100 g
Rear beam	100 g x 10 g
Stainless steel plate platform	5.8 x 5.8 inches / 14.73 x 14.73 cm

## Accessories

	Ohaus Number
Security device for loop .....	76288-01
Security device for Kensington slot .....	470004-010
Replacement attachment weight set .....	707-00
Replacement rod & clamp assembly .....	183-00

## OHAUS CORPORATION

\* 19A Chapin Road  
P.O. Box 2033  
Pine Brook, NJ 07058 USA

Tel: 800.672.7722  
973.377.9000  
Fax: 973.944.7177

www.ohaus.com

With offices throughout Europe,  
Asia, and Latin America

\* ISO 9001:2000  
Registered Quality Management System

80774440 © Copyright Ohaus Corporation

